

UNICSKIN

YOUR EFFECTIVE SOLUTION

UNICLED NECK & CHIN 5.0

FIRMER NECK · DEFINED JAWLINE · EVEN TONE

CLINICAL TECHNOLOGY · PROFESSIONAL RESULTS AT HOME



**LONGEVITY SCIENCE
FOR THE NECK, JAWLINE AND DÉCOLLETÉ**



LED light
GREEN

SLIMS DOWN
DOUBLE CHIN




LED light
RED

SMOOTHS
WRINKLES



LED light
DEEP RED

EVEN'S OUT
DARK SPOTS

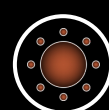

LIGHT
NIR

FIRMS THE NECK &
DÉCOLLETAGE

USER MANUAL



LET'S GET STARTED

Congratulations! You have taken the first step towards a firmer, smoother, and visibly rejuvenated neck and chin contour, with UNICSKIN, global experts in medical beauty technology and at-home skincare.

Before you start using this device, **carefully read all instructions and precautions** and use it only for the purpose described in this manual.

1 DEVICE, 2 TECHNOLOGIES, 2 MODES, 6 BENEFITS

UNICLED NECK & CHIN 5.0 represents a new generation of smart skincare for the neck and chin. Its technology combines **red LED light, green LED light, and NIR light** in a wireless, easy-to-use device designed to **integrate seamlessly into your daily wellness and longevity routine**.

In the morning, it helps to revitalize the skin and promote a brighter and more even appearance. At night, it helps to smooth wrinkles and expression lines, as well as visibly improving skin firmness and elasticity.

MAIN BENEFITS



SLIMS DOWN
DOUBLE CHIN



REAFFIRMS
THE NECK



EVENS OUT SKIN
TONE



REDEFINES THE
FACIAL OVAL



SMOOTHS
WRINKLES



REDUCES DARK
SPOTS

THE SCIENCE BEHIND A YOUNGER-LOOKING NECK

This device combines green, red and intense red LED light with NIR light to treat the **neck, jaw and décolleté in a comprehensive way**.

Each technology works in a complementary way to help treat the main visible signs of aging in the area: double chin, sagging, wrinkles, lines, uneven tone and dark spots.

Each technology has a specific function:



GREEN LED LIGHT (520nm)

It acts primarily on the epidermis, where melanin is produced. It helps regulate melanin synthesis to **reduce the appearance of dark spots and hyperpigmentation**, providing a more even and radiant complexion. It is part of the **WRINKLES + DARK SPOTS Mode** for the neck and décolleté area.

UNISKIN

YOUR EFFECTIVE SOLUTION

RED LED LIGHT (630nm)

It works in the dermis, where it **stimulates fibroblast activity and the natural synthesis of collagen and elastin**, helping to reduce collagenase (the enzyme responsible for aging). It provides **greater firmness, elasticity, and smoothness** to the neck and jawline. It is part of the **SLIM** and **WRINKLES + DARK SPOTS Modes**.

INTENSE RED LED LIGHT (660nm)

It penetrates deeper into the dermis to **enhance the effects of red light**, helping to smooth the appearance of deep wrinkles and firm the skin's structure. It is part of the **WRINKLES + DARK SPOTS Mode** for the neck and décolleté area.

NEAR INFRARED LIGHT (850nm)

It has the longest wavelength and reaches the hypodermis, the deepest layer of the skin. It helps **activate cellular energy and deep regeneration**, and stimulates the metabolism of adipocytes (the cells that store fat), contributing to firming and redefining the contour of the chin and jawline. It is part of the **SLIM** and **WRINKLES + DARK SPOTS Modes**.

In addition, **LED + NIR light** enhances microcirculation and cellular oxygenation, improving the absorption and performance of dermocosmetics applied afterwards.

TWO MODES. TWO SPECIFIC TREATMENTS.



SLIM MODE

Chin Reduction + Facial Oval Definition

Uses **Red LED** and **NIR** light for:

- Stimulate the metabolism of adipocytes in the chin area.
- Reaffirm and lift the jawline.
- Redefine facial contours and firm.

Treatment areas: Double chin and jawline.

WRINKLES/DARK SPOTS MODE

Anti-aging + Depigmenting

Combines **Red**, **Intense Red**, **Green** and **NIR** LED lights to:

- Reduce wrinkles and expression lines.
- Add luminosity and unify skin tone.
- Minimize dark spots.

Treatment areas: Neck and décolleté.



UNICLED NECK & CHIN 5.0 USER GUIDE

SLIM MODE

Press and hold the button for 1 second to turn on SLIM Mode (double chin and jaw).

WRINKLES + DARK SPOTS MODE

Press and hold the button for 1 second to turn on WRINKLES + DARK SPOTS Mode (neck and décolletage).

MEDICAL GRADE ANTIBACTERIAL SILICONE

Soft, flexible, and antibacterial. It adapts to every area without leaving marks and can be safely shared.

USB CHARGING PORT


100% wireless. It recharges via USB and a single charge allows for multiple sessions.

ERGONOMIC AND LIGHTWEIGHT DESIGN

It adapts perfectly to the contour of the neck and chin, so you can use it while performing other tasks.


Thanks to its ergonomic design, the UNICLED NECK & CHIN 5.0 directs the LED/NIR light towards the neck, jawline, and décolleté. **Its shape adapts to the area's curvature for even treatment, respecting the delicate skin.**

UNICLED NECK & CHIN 5.0 USAGE INSTRUCTIONS

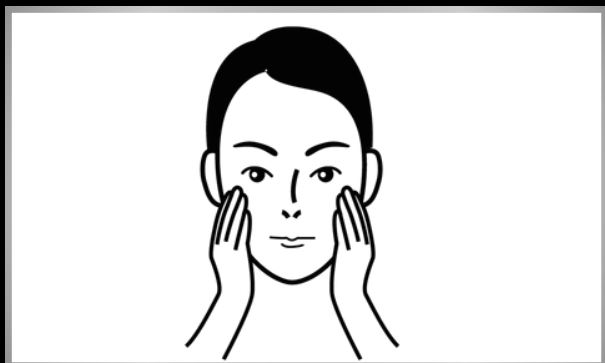


STEP 1

Cleanse and dry the skin of your face, neck, and décolleté before starting. The treatment should always be applied to clean skin, before your beauty routine, so that the light can work at its most effective.

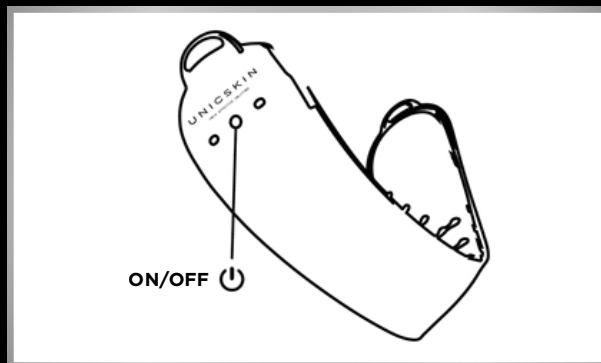

STEP 2

Turn on the device by pressing and holding the button for 1 second. Make sure it is fully charged: you will know when the indicator lights remain solid green.

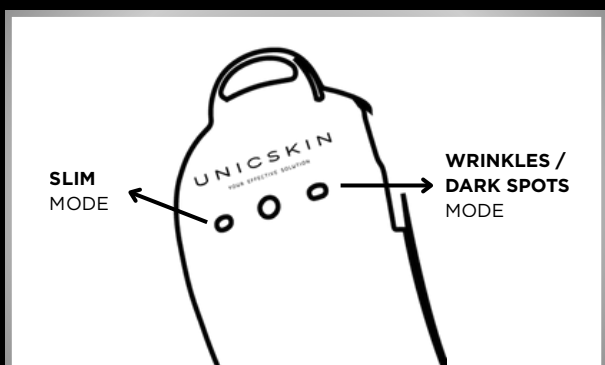

STEP 3

Select the mode according to the area you want to treat: WRINKLES + DARK SPOTS for the neck and décolleté, or SLIM for the double chin and jawline.



STEP 4

Place the device on the chosen area and leave it on for 20 minutes. It will automatically turn off at the end of the session. Then, continue with your routine.

For optimal results, we recommend using UNICLED NECK & CHIN 5.0 for **20 minutes** per session, a **minimum of 3 times per week for 4 weeks**. To maintain results, continue 1 or 2 times per week.

Complete your routine: in **SLIM Mode**, with the UNICLIFTING+ firming serum; in **WRINKLES + DARK SPOTS Mode**, with the UNICWRINKLE 360 serum and the UNICD-PIGMENT XTREME COMPLEX depigmenting cream.

UNICLED NECK & CHIN 5.0 ACCESSORIES


UNICLED NECK & CHIN 5.0 comes in a compact and elegant case that includes everything you need to start treatment comfortably, safely and efficiently.



UNICLED NECK
& CHIN 5.0
device



Travel case

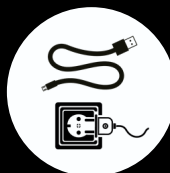

USB-C charging
cable



Express user
manual

LOAD

The UNICLED NECK & CHIN 5.0 includes its USB Type-C charging cable. A full charge (approx. 3-4 h) allows up to 6 sessions.

HOW TO CARE FOR THE DEVICE

After each use, gently clean the surface of the UNICLED NECK & CHIN 5.0, especially the areas that have been in contact with your skin. Use a soft, slightly damp cloth or a water-based, alcohol- and perfume-free cleaning wipe, and allow the device to dry completely before storing it.

Always store the device in its neoprene case to protect it from dust, scratches and possible deformation, and to keep it in optimal condition between uses.

TECHNICAL SPECIFICATIONS

Power	5 W
Voltage	5 V
Irradiance SLIM	15 mW/cm ²
Irradiance W+DS	25 mW/cm ²
Wavelengths	520 / 630 / 660 / 850 nm
N° of LED lights	112
Material	Antibacterial medical- grade silicone
Conservation	Cool place, max. 42 °C



PRECAUTIONS

- Do not look directly at the lights or place the device over your eyes.
- The device is not waterproof. Do not submerge, wash under running water, or get the charging port wet.
- Do not disassemble, tamper with, or attempt to repair the device.
- Not recommended for children under 12, pregnant or breastfeeding women, people sensitive to light, or those with medical conditions. Consult your doctor if you have a thyroid condition, a pacemaker, metal implants, or active cancer.
- Keep out of reach of children.

FREQUENTLY ASKED QUESTIONS

How does LED + NIR technology work on the neck and double chin?

In SLIM Mode, red light (630nm) and NIR (850nm) help promote microcirculation, tissue metabolism, and firmness. Since a double chin is related to submental fat, the chin appears slimmer and the contour more defined over time. In WRINKLES + DARK SPOTS Mode, green light (520nm), red light (630nm), intense red light (660nm), and NIR (850nm) help promote collagen production, improve texture, and soften the appearance of wrinkles and dark spots on the neck and décolleté.

What's new in version 5.0?

The WRINKLES + DARK SPOTS mode now incorporates a new green LED light (520nm) to reinforce its action on spots and tone uniformity.

How long does the session last and how often should I use it?

Each session lasts 20 minutes and the device turns off automatically. Use it 3 times a week for 4 weeks, then 1-2 times a week to maintain results.

Does the treatment hurt?

No. LED + NIR phototherapy is non-invasive and painless. It has CE and FDA certifications.

Can I treat my neck and double chin on the same day?

Yes, following the recommended 20-minute session for each mode.

Is it safe in the thyroid area?

It uses non-invasive LED + NIR phototherapy. As this is a sensitive area, if you have a thyroid condition or are under medical treatment, consult your doctor before using it.

Can I use it during pregnancy or breastfeeding?

Although LED and NIR light are non-invasive and well-tolerated, their use has not been specifically evaluated during pregnancy or breastfeeding, so we recommend consulting your doctor.

CUSTOMER SERVICE

If you have any questions or comments you would like to share with us, please write to us at the following email address: customerservices@unicskin.com.

DISCLAIMER

This product is not intended to diagnose, cure, or prevent any disease. UNICSKIN S.L.

This product does not guarantee its ability to treat any physical, mental, or skin condition. For any health concerns, consult a qualified medical professional.

WARRANTY

Except in cases of manufacturing defects, returns are not accepted for products purchased through our website (www.unicskin.com) once they have been used and opened. If the product was purchased outside of our website, please consult the warranty policy offered by the retailer, which in no case supersedes this warranty policy.

If the packaging has not been opened, UNICSKIN will accept returns within 14 days of the purchase date. The customer is responsible for all return shipping costs and risks. The device must be returned in its original packaging, along with the original proof of purchase. Please contact Customer Service at customerservices@unicskin.com or by phone at +34 912 797 377 for information on the return address.

In the event of a manufacturing defect and during the first three years from the date of purchase, the return of any of our Techno-Beauty products may result in a replacement with the same product (in which case the return costs will be borne by UNICSKIN), a refund, or a price reduction if the defect detected is of little importance or the value of the item at the time of the claim does not justify its replacement due to the time elapsed.

In any case, it will be necessary to keep and present the receipt or invoice for the purchased product, clearly showing the purchase date. Furthermore, UNICSKIN reserves the right to carry out the necessary checks on the purchased product, requesting photographs or evidence of the alleged nonconformity or the return of the defective product. Initially, the customer will be responsible for the shipping costs, which will be refunded once UNICSKIN has verified that the detected fault corresponds to a manufacturing defect.

Products damaged due to storage conditions, maintenance, or any other type of handling by the customer are excluded from this warranty policy. The product complies with the legal and regulatory standards applicable in the European Community at the time of manufacture and distribution.

UNICSKIN

YOUR EFFECTIVE SOLUTION

GLOBAL EXPERTS IN TECHNO BEAUTY & HOME BEAUTY
EFFECTIVE RESULTS PROVEN BY SCIENCE
SUSTAINABILITY & FORMULAS +92% OF NATURAL ORIGIN
INNOVATIVE LUXURY WITHIN YOUR REACH
