

UNICSKIN

YOUR EFFECTIVE SOLUTION



UNICBLUE

UNICBLACK

UNICGREEN

UNICPINK

UNIC**THERMO-SONIC** 2.0 **FUN COLLECTION**

MULTI-ACTION FACIAL BRUSH
THE MOST VERSATILE ANTI-AGING BEAUTY DEVICE

FOR ANYONE. ANYWHERE. ANYTIME

UNIC**DEEP**
CLEANSE

UNIC**FACE**
YOGA

UNIC**EYES**
& LIPS



UNICKSKIN

YOUR EFFECTIVE SOLUTION

1 DEVICE, 2 TECHNOLOGIES, 3 ACTIONS



HEAT TECHNOLOGY

Heat treatment at 42° to reduce the signs of aging, for visibly younger looking skin.

SONIC TECHNOLOGY

10,000 pulses per minute for powerful lymphatic drainage that decongests, for visible more radiant skin.

TECHNOLOGY THAT **MULTIPLIES x3**
THE EFFECTS OF AN ACTIVE DERMOCOSMETICS

APPROVED BY EXPERTS AS A SAFE AND EFFECTIVE SKIN CARE DEVICE

UNICSKIN


YOUR EFFECTIVE SOLUTION



1
**UNICDEEP
CLEANSE**



2
**UNICFACE
YOGA**

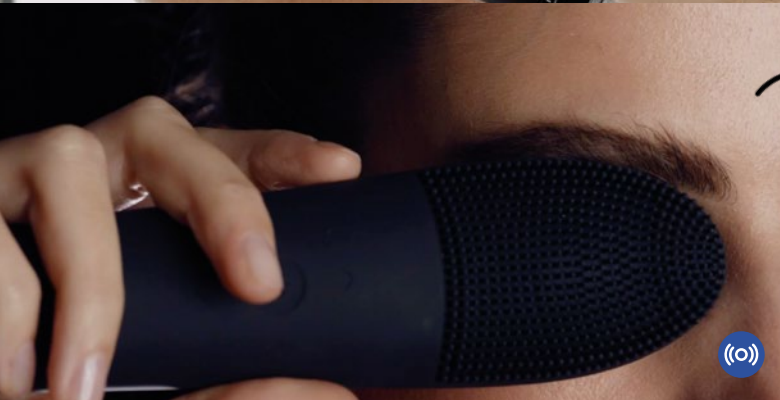

3
**UNICEYES
& LIPS**



1 DEVICE, 2 TECHNOLOGIES, 3 ACTIONS

UNICSKIN

YOUR EFFECTIVE SOLUTION

1 DEVICE, 2 TECHNOLOGIES, 3 ACTIONS



UNICDEEP
CLEANSE

Deep cleanse, make-up removal and anti-pollution action

 **42 C° Anti-aging Heat Technology + Sonic Technology with 10,000 pulses per minute.**

With a smooth exfoliation suitable for all skin types (including the most sensitive), flushing out pores and eradicating acne-causing bacteria.



UNICFACE
YOGA

Anti-aging Treatment (UNICFACE YOGA)

 **42 C° Anti-aging Heat Technology + Sonic Technology with 10,000 pulses per minute.**

Significantly improves skin elasticity and smooths out wrinkles with our facial gym, leaving you with softer, velvety skin.



UNICEYES
& LIPS

Eye and lip contour treatment

 **Intermittent Sonic Technology (Lymphatic Drainage).**

Treatment for dark circles, bags and wrinkles through a powerful lymphatic drainage. Lip contour treatment smoothing wrinkles and lines around the mouth. A massage with Technology that optimizes the results of manual massages.

YOUR BEST PARTNER FOR A YOUNG AND RADIANT SKIN!

UNICSKIN

YOUR EFFECTIVE SOLUTION

MODE OF USE

1 UNICDEEP CLEANSE


DEEP CLEANSE, MAKE-UP REMOVAL AND ANTI-POLLUTION ACTION (100% Waterproof)

 **42 C° Anti-aging Heat Technology + Sonic Technology with 10,000 pulses per minute.**

Recommended to use the brush with UNICSKIN product: UNICDETOX&GO All-in-One cleanser and toner foam.




Moisturize your face

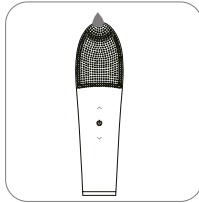


Add the foam on the brush head



Turn on the Heat and Sonic Technology by pressing once

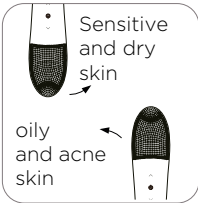

Choose the best option for your skin needs

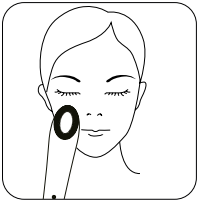

Cleanse the face in circular movements



Reach the most difficult areas



 **30 SEG.**



Clinically proven results: more than **93%** of the volunteers saw their pores cleaner and their skin clearer



UNICSKIN

YOUR EFFECTIVE SOLUTION

2 UNICFACE YOGA


ANTI-AGING TREATMENT FOR MATURE SKIN

 **42 C° Anti-aging Heat Technology + Sonic Technology with 10,000 pulses per minute.**

Recommended to use the brush with UNICSKIN's products: UNIC**Serums** (UNIC**HYDRO-ACTI-VE**, UNIC**LIFTING+**, UNIC**LUMINOSITY 3.0** y UNIC**WRINKLE-STOP**).




Apply the serum
on your face...

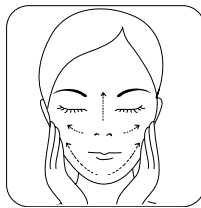

... with upward
movements




Work the forehead
area

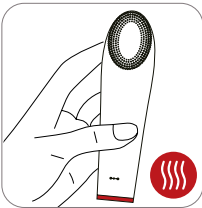

Turn on the Heat and
Sonic Technology by
pressing once



Massage the area
to be treated with
the brush

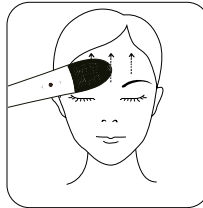

and relax the
muscles



 **3 MIN.**

“BIKRAM YOGA” FOR YOUR FACE: facial yoga massage with **Heat Technology at 42°** equivalent to an anti aging spa treatment in only 3 minutes



UNICFACE YOGA

THE TECHNOLOGICAL REVOLUTION OF FACIAL GYMNASTICS

The **42°C Heat Technology** of our **UNICTHERMO-SONIC 2.0** facial brush **brings an anti-aging spa treatment to every home in just 3 minutes**. **UNICSKIN's unique UNICFACE YOGA exercises multiply the effects of manual facial gymnastics, achieving immediate results.**

UNICSKIN protocol: Every day, practice 5 facial gymnastics exercises that will help you tone your face, keep it tight, juicy and prevent skin sagging. Include this routine along with your abdominal or squats series.



UNICFACE YOGA LIPS

Use the back side of the brush in gentle circular movements around the lips with mode 1 (Heat Technology and Sonic Technology) together with the serum suitable for your skin needs. At the same time, run your tongue around the inside of the lips in the same direction 3 times, and repeat in the opposite direction. Swallow. **This helps to smooth nasolabial folds and fine lines around the mouth, lifting and strengthening the muscles.**




UNICFACE YOGA RADIANT YOUNG SKIN

Gently press the short silicone bristles to your skin for 2" and glide up the brush with the back side in mode 1 (Heat Technology and Sonic Technology) to your cheeks, forehead, nasolabial fold, and chin for 60" with your regular UNICSKIN serum. **This exercise improves the blood circulation and stimulates the production of volume agents. It is a perfect technique to make your skin glow in the morning.**



UNICFACE YOGA EYE LIFTING

Gently press the back side of the brush in mode 1 (Heat Technology and Sonic Technology) on the outer tail of the eyebrows. Open the eyes wide and stretch the eyebrows with the brush. Hold the position for 10" and repeat 2 more times. **This exercise reduces the lines between the eyebrows, relieves tension and firms the forehead skin. With UNICSKIN Technology, you will multiply the effects of a manual facial yoga exercise.**

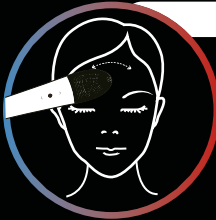

UNICFACE YOGA FOREHEAD WRINKLES

Gently press the brush to the center of the forehead, with the brush's back side, and glide the device across your forehead sides, for 3" each side, in mode 1 (Heat Technology and Sonic Technology) together with the serum of your choice. Repeat 3 times. **This gentle massage technique improves the circulation, increases the collagen and the elastin production, and releases the tension that can cause deep wrinkles.**



UNICFACE YOGA CHEEKBONES ELEVATION

Lightly press the back side of the brush in mode 1 (Heat Technology and Sonic Technology) under the cheekbone and gently glide the device while pressing and lifting the muscle. Hold the lift for 5" and repeat the movement. Follow the path of the bone supporting the cheeks until you reach the temple. Repeat the same exercise with the opposite cheekbone.

UNICSKIN

YOUR EFFECTIVE SOLUTION

3 UNICEYES & LIPS

EYE AND LIP CONTOUR TREATMENT

 Intermittent Sonic Technology (Lymphatic Drainage).


Recommended to use the brush with UNICSKIN's products: UNICEYES TRIPLE ACTION Eye and Lip Contour + collagen patches UNICCOLLAGEN FLASH EYE MASK or UNICCOLLAGEN VOLUME LIP MASK.



Apply eye and lip contour



Turn on the Intermittent Sonic Technology by pressing twice

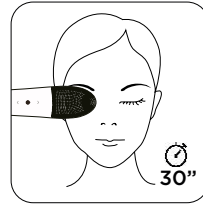

Hold and gently move in circles the brush over each eye



Massage the lip contour



Add our collagen eye patch



Add our collagen lip patch



 60 SEG.




A massage with Technology that optimizes the results of manual massages.

UNICSKIN

YOUR EFFECTIVE SOLUTION

UNIC THERMO-SONIC 2.0 FUN COLLECTION OVERVIEW

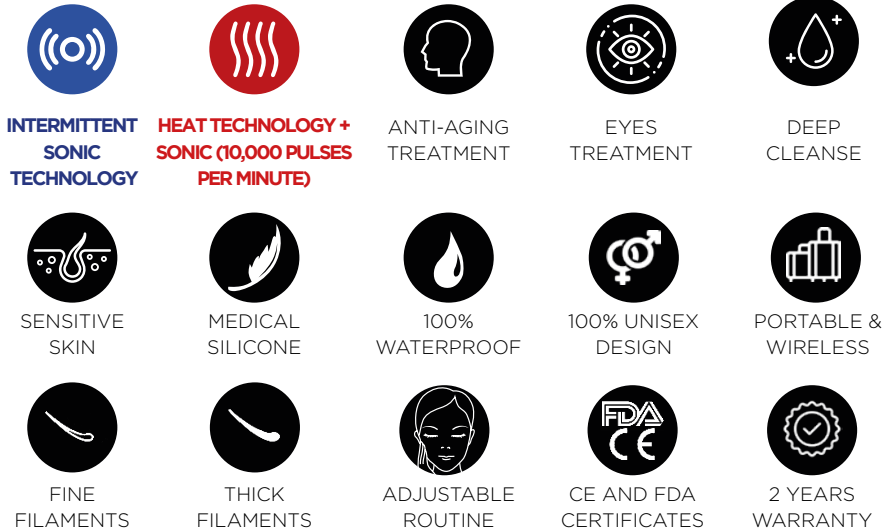

THE TECHNOLOGICAL REVOLUTION OF FACIAL GYMNASTICS

YOUR DAILY DOSE OF FUN!



READY FOR A NEW SKINCARE ERA?





YOUR EFFECTIVE SOLUTION



MULTI-ACTION FACIAL BRUSH
AVAILABLE IN 4 DIFFERENT COLORS,
DEPENDING ON YOUR OWN PERSONALITY,
ALL WITH AN ECO-FRIENDLY PACKAGING
LOVING OUR WORLD.
LOVING YOU.

UNICSKIN

YOUR EFFECTIVE SOLUTION

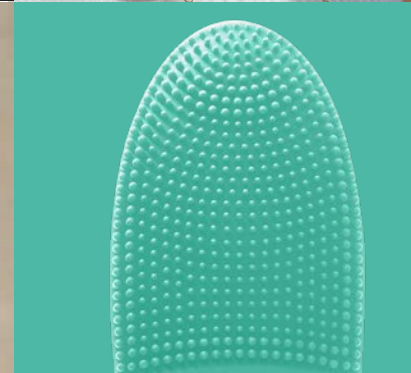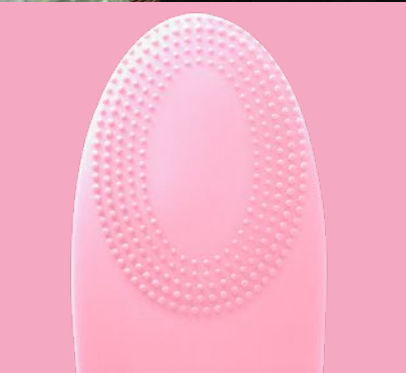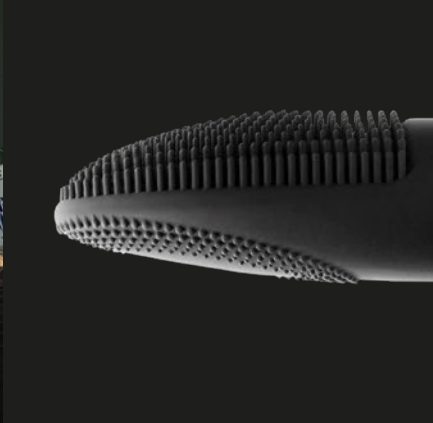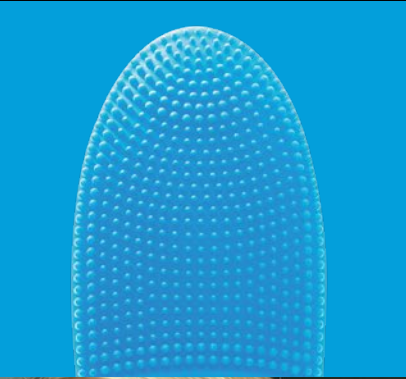

www.unicskin.com
customerservices@unicskin.com

#UNICSKINERS

